



Pipette fillers, Easypet®



eppendorf

The device responds exactly to finger pressure for slow or rapid aspiration, setting the meniscus or emptying the pipette.

- Active pressure compensation prevents corrosion caused by aspirated vapour
- High aspirating speed: 25mL in 5s
- One finger control for the flow out and blow out functions
- Cordless work for up to eight hours, rechargeable
- Autoclavable pipette adapters and filter holders
- Space saving wall holder


Precise, safe and fatigue free series pipetting, using graduated pipettes and one mark pipettes, is possible with this light, hand held device.

Suitable for pipettes of 0.1mL to 100mL. Non-slip, even with wet pipettes. Smooth setting of the pumping efficiency. Fine control of the meniscus movement. With two 0.45µm membrane filters, recharger and wall holder. Autoclavable adapters and filter holders.

Catalogue No	Legacy No	Alt. No	Description
10759248	-	4420000036	Easypet® with transformer, wall holder and two membrane filters/0.45µm

Spare parts

Catalogue No	Legacy No	Alt. No	Description	Pack qty
10296571	-	4420803007	Membrane filter, 0.45µm, sterile	5
11604412	-	4421601009	Membrane filter, 0.45µm, PTFE, sterile	5
10328611	-	4420810003	NiCd batteries	3

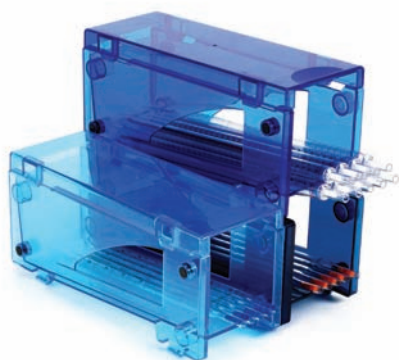

Pipette rack, stackable, Fisherbrand®



Fisherbrand®

Ideal for storing and organising serological pipettes.

- Angled shelves store 1.0mL to 50mL pipettes within easy reach
- Click together to horizontally or lock in place with top connectors
- Protective lid to prevent dust build up
- Magnets on both sides of the unit connect securely and safely to metal surfaces
- Material: Acrylonitrile butadiene styrene (ABS)

Catalogue No	Legacy No	Description	Dimensions [l x w x h], mm
11958004	-	Supplied as 3 compartments with 3 lids	286 x 301 x 127.5



Pipette support rack

scienceware
BET-ART PRODUCTS

For orderly, space saving pipette storage, these sturdy racks hold 50 pipettes in a vertical 5 x 10 arrangement. Top and middle decks have 16mm dia. holes; bottom deck holes are tapered from 16mm to 6mm dia. Pipette tips will not contact table surface. The end plates have handles for easy carrying. Autoclavable at 121°C.

Catalogue No	Legacy No	Alt. No	Width, mm	Depth, mm	Height, mm
10223143	-	F18953-0000	213	114	222